

# Rapid Detection of Polymyxin Resistance in *Enterobacteriaceae*

Patrice Nordmann, Aurélie Jayol, Laurent Poirel

For identification of polymyxin resistance in *Enterobacteriaceae*, we developed a rapid test that detects glucose metabolism associated with bacterial growth in the presence of a defined concentration of colistin or polymyxin B. Formation of acid metabolites is evidenced by a color change (orange to yellow) of a pH indicator (red phenol). To evaluate the test, we used bacterial colonies of 135 isolates expressing various mechanisms of colistin resistance (intrinsic, chromosomally encoded, and plasmid-mediated MCR-1) and 65 colistin-susceptible isolates. Sensitivity and specificity were 99.3% and 95.4%, respectively, compared with the standard broth microdilution method. This new test is inexpensive, easy to perform, sensitive, specific, and can be completed in <2 hours. It could be useful in countries facing endemic spread of carbapenemase producers and for which polymyxins are last-resort drugs.

Among the most clinically significant multidrug-resistant bacteria are carbapenemase-producing *Enterobacteriaceae*. Because these bacteria usually remain susceptible to polymyxins, an old class of antimicrobial drugs almost abandoned in the 1970s because of their potential toxicity, interest in polymyxins (colistin and polymyxin B) has been renewed worldwide (1,2). However, the increasing use of colistin explains why acquired colistin resistance may now be added to the carbapenem resistance trait in *Enterobacteriaceae* (3).

The standard reference technique for determining susceptibility to polymyxins is broth microdilution, which requires fastidious attention and a long time (24 h) to perform (4). Other techniques for determining susceptibility to polymyxins (disk diffusion and Etest) have been proposed and also provide results in 18–24 h. Because of poor diffusion of polymyxin molecules in agar, rates of false susceptibility are high (up to 32%) (4,5).

Acquired resistance to colistin in *Enterobacteriaceae* results mostly from modification of lipopolysaccharide (6). Addition of phosphoethanolamine, 4-amino-L-arabinose cationic groups, or both to lipopolysaccharide decreases polymyxin binding to the bacterial outer membrane. Addition of these groups may be associated with chromosome-encoded mechanisms (mutations in PmrAB or PhoPQ

two-component systems or alterations of the *mgrB* gene) (6). A recent report revealed that addition of phosphoethanolamine may also be plasmid mediated through the *mcr-I* gene, which confers the first known plasmid-mediated resistance to colistin in isolates from humans and animals (7). More recently, the *mcr-I* gene was identified in several plasmid backbones, mostly in *Escherichia coli* (8–10). There is therefore a need for a test that enables rapid detection of polymyxin resistance in *Enterobacteriaceae* and that may contribute to its containment.

We developed a test (the rapid polymyxin NP [Nordmann/Poirel] test) that detects bacterial growth in the presence of a defined concentration of a polymyxin. Bacterial growth detection (or absence) is based on carbohydrate metabolism (11). Acid formation associated with carbohydrate metabolism in *Enterobacteriaceae* can be observed through the color change of a pH indicator. This test is rapid (<2 h) and easy to perform.

## Materials and Methods

### Isolate Collection

To evaluate the performance of the rapid polymyxin NP test, we used 200 isolates collected from clinical samples worldwide. This collection included 135 *Enterobacteriaceae* isolates resistant to polymyxin: 5 isolates of intrinsically polymyxin-resistant species (*Morganella morganii*, *Proteus mirabilis*, *Proteus vulgaris*, *Providencia stuartii*, and *Serratia marcescens*) and 130 isolates of various enterobacterial species (*Klebsiella* spp., *E. coli*, *Enterobacter* spp., and *Hafnia alvei*) with acquired resistance to polymyxins (online Technical Appendix, <http://wwwnc.cdc.gov/EID/article/22/6/15-1840-Techapp1.pdf>), including a previously reported heteroresistant *Klebsiella pneumoniae* isolate for which MIC for colistin was high (12). The other 65 enterobacterial isolates belonged to various species and were susceptible to polymyxins (online Technical Appendix).

### MIC Determination

To determine MICs for polymyxins, we used the broth microdilution method in cation-adjusted Mueller-Hinton broth (MHB-CA, reference 69444; Bio-Rad, Marnes-La-Coquette, France) as recommended by Clinical Laboratory

Author affiliation: University of Fribourg, Fribourg, Switzerland

DOI: <http://dx.doi.org/10.3201/eid2206.151840>

Standard Institute (CLSI) guidelines (13–15). We considered this method to be the standard for comparison with the rapid polymyxin NP results. Polymyxin B and colistin sulfate (Sigma-Aldrich, St. Louis, MO, USA) were tested over a range of dilutions (0.12–128 µg/mL). All experiments were repeated in triplicate in separate experiments. As recommended by CLSI, microdilution was performed without addition of Tween 80 (15), and *E. coli* ATCC 25922 was used as a control strain.

Because no breakpoint is available for polymyxins for *Enterobacteriaceae* according to CLSI guidelines (14), we used the breakpoints of the European Committee on Antimicrobial Susceptibility Testing (EUCAST) for reference (16). Enterobacterial isolates with colistin or polymyxin B MICs  $\leq 2$  µg/mL were categorized as susceptible; those with MICs  $>2$  µg/mL were categorized as resistant.

### PCR Amplification and Sequencing

We recovered the chromosomal DNA of the isolates by using a commercially available kit (QIAquick; QIAGEN, Courtaboeuf, France) according to the manufacturer's instructions. We sequenced the *pmrA*, *pmrB*, *phoP*, *phoQ*, and *mgrB* genes possibly involved in colistin resistance in *K. pneumoniae* and *K. oxytoca*, as described previously (12,17–19). We performed PCR amplification for detection of the plasmid-mediated *mcr-1* gene as described (7).

### Isolate Genotyping by Pulsed-Field Gel Electrophoresis

We assessed the genetic relatedness of the colistin-resistant isolates with identical molecular mechanisms of colistin resistance. We used pulsed-field gel electrophoresis with *Xba*I-digested genomic DNA as described previously (20).

### Rapid Polymyxin NP Test

#### Reagents and Solutions

The rapid polymyxin NP test uses 2 reagents and solutions: stock solutions of polymyxins and rapid polymyxin NP solution. Each is described below.

For stock solutions of polymyxins, colistin sulfate and polymyxin B powders (Sigma Aldrich) were diluted into MHB-CA medium in glass tubes to obtain a concentration of 0.2 mg/mL. These powders can be stored at 4°C before use, and the diluted polymyxin solutions can be kept at –20°C for 1 year. Of note, polymyxin-containing batches from commercial origin can be used, but the colistimethate sulfate powder, a therapeutic prodrug of colistin, cannot be used.

To prepare 250 mL of the rapid polymyxin NP solution, we mixed the culture medium and the pH indicator in a glass bottle as follows: 6.25 g of MHB-CA powder, 0.0125 g of phenol red (Sigma Aldrich), and 225 mL of distilled water. The pH of the solution was adjusted to 6.7

by adding drops of 1 mol/L HCL. This solution was then autoclaved at 121°C for 15 min. After cooling the solution to room temperature, we added 25 mL of 10% anhydrous D(+)-glucose (Roth, Karlsruhe, Germany) sterilized by filtration. The final concentrations in the rapid polymyxin NP solution were consequently 2.5% of MHB-CA powder, 0.005% of phenol red indicator, and 1% of D(+)-glucose. This rapid polymyxin NP solution can be kept at 4°C for 1 week or at –20°C for 1 year. This solution must be pre-warmed at 37°C before use to prevent growth delay and therefore a delayed color change.

Just before performing the experiment, we added colistin to the rapid polymyxin NP solution and mixed it into sterile glass tubes to obtain a rapid polymyxin NP solution containing a final colistin concentration of 5 µg/mL. For example, we added 25 mL of colistin stock solution at 0.2 mg/mL to 1 mL of rapid polymyxin NP test solution for the testing of 1 isolate and respective negative and positive controls.

#### Bacterial Inoculum Preparation

We prepared a standardized enterobacterial inoculum by using freshly obtained (overnight) bacterial colonies grown on Luria-Bertani or Mueller-Hinton plates. We resuspended the bacterial colonies into 10 mL of sterile NaCl (0.85%) to obtain a 3.0–3.5 McFarland optical density ( $\approx 10^9$  CFU/mL), which corresponds to an  $\approx 10$ -mL full loop of bacterial colonies diluted in 10 mL of NaCl. A bacterial suspension was prepared for each isolate to be tested and for the colistin-susceptible and -resistant isolates used as controls (isolates FR-01 and FR-136; respectively; online Technical Appendix). As recommended by the EUCAST guidelines for susceptibility testing, we used the bacterial suspensions within 15 min of preparation and for no longer than 60 min after preparation (16).

#### Tray Inoculation

We performed testing in a 96-well polystyrene microtest plate (round base, with lid, sterile, reference 82.1582.001; Sarstedt, Nümbrecht, Germany). For each isolate, bacterial suspension was inoculated in parallel into 2 wells, with and without colistin. The following steps of the rapid polymyxin NP test were then performed (Figure):

Step 1: 150 µL of colistin-free solution was transferred to wells A1–A4.

Step 2: 150 µL of the rapid polymyxin NP solution containing colistin was transferred to wells B1–B4.

Step 3: 50 µL of NaCl 0.85% was added to wells A1 and B1.

Step 4: 50 µL of the colistin-susceptible isolate suspension used as negative control was added to wells A2 and B2.

Step 5: 50 mL of the colistin-resistant isolate suspension used as positive control was added to wells A3 and B3.

Step 6: 50 mL of the tested isolate suspension was added to wells A4 and B4.

We mixed the bacterial suspension with the reactive medium by pipetting up and down. The final concentration of bacteria was  $\approx 10^8$  CFU/mL in each well, and the final concentration of colistin sulfate was 3.75  $\mu$ g/mL.

#### Tray Incubation

We incubated the inoculated tray for up to 4 h at  $35 \pm 2^\circ\text{C}$  in ambient air, without being sealed and without agitation. We did not seal the tray because oxygen is required for carbohydrate metabolism.

#### Tray Reading

We visually inspected the tray (checked for no spontaneous color change) after 10 min and then every hour for 4 h. We considered the test result positive (polymyxin resistance) if the polymyxin-resistant isolate grew in presence of colistin and negative (polymyxin susceptibility) if the polymyxin-susceptible isolate did not grow in presence of colistin. We considered the test result interpretable if the following 4 conditions were met: 1) both wells with 0.85% NaCl without bacterial suspension (wells A1 and B1) remained orange (absence of medium contamination); 2) the wells with bacterial suspension and colistin-free (wells A2–A4) turned from orange to yellow, confirming the metabolism of glucose by the isolates; 3) the wells with the colistin-susceptible reference bacterial suspension (negative control)

gave negative results (wells A2 and B2); and 4) the wells with the colistin-resistant reference bacterial suspension (positive control) gave positive results (wells A3 and B3). The test result was positive when the well containing colistin (well B3) and the isolate to be tested turned from orange to yellow, giving exactly the same color as the well without colistin (well A3), indicating glucose metabolism and growth in presence of colistin (i.e., colistin resistance) (Figure). The test result was negative when the well containing colistin (well B2) with the isolate to be tested remained orange (unchanged color) (Figure) or was more clear than the wells with 0.85% NaCl but not exactly the same color as the well without colistin (not shown). Results were interpreted by 2 technicians who did not know which isolates were colistin resistant and colistin susceptible.

#### Other Experimental Conditions Tested

##### Polymyxin B Instead of Colistin

We evaluated the possibility of adapting the test to susceptibility testing of polymyxin B in countries where polymyxin B is prescribed. To do so, we performed the rapid polymyxin NP test with 20 colistin/polymyxin B-susceptible isolates and 20 colistin/polymyxin B-resistant isolates with polymyxin B at the same concentrations of colistin sulfate.

##### Incubation Conditions


To determine effects of the incubation atmosphere on time to result, we incubated the tray with 20 colistin-susceptible isolates and 20 colistin-resistant isolates, in parallel, under 2 conditions: ambient air and atmosphere with 5% CO<sub>2</sub>. We also incubated the tray for 20 colistin-susceptible isolates and 20 colistin-resistant isolates in parallel with and without agitation.

##### Culture Media

To determine the potential effects of culture medium on the test results, we performed the test with 20 colistin-susceptible isolates and 20 colistin-resistant isolates cultured overnight on different agar plates. The following culture media were tested: 1) nonselective culture medium such as Columbia agar + 5% sheep blood (bioMérieux, La-Balme-Les-Grottes, France); 2) chocolate agar + PolyVitex (bioMérieux); 3) nonselective chromogenic medium UriSelect 4 (Bio-Rad); 4) Eosin methylene blue agar (Sigma Aldrich); 5) Drigalski agar (Bio-Rad); 6) MacConkey agar (VWR BDH Prolabo, Leuven, Belgium); and 7) bromocresol purple (bioMérieux).

#### Results

Of the 200 enterobacterial isolates tested to evaluate the performance of the rapid polymyxin NP test (online Technical Appendix), 5 isolates belonged to bacterial species that are

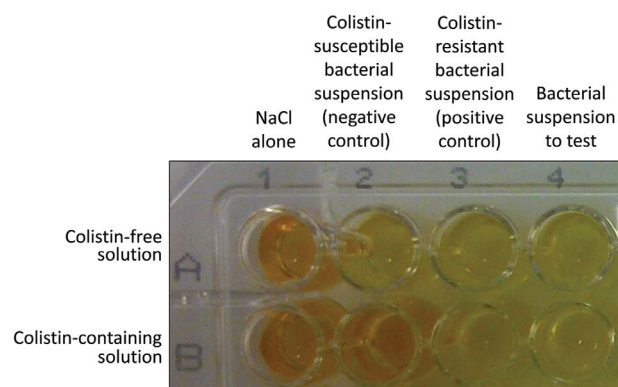

**Figure.** Representative results of the rapid polymyxin NP [Nordmann/Poirel] test. Noninoculated wells are shown as controls (first column). The rapid polymyxin NP test was performed with a reference colistin-susceptible isolate (second column) and with a reference colistin-resistant isolate (third column) in a reaction medium without (upper row) and with (lower row) colistin. The tested isolate grew in the presence (and absence) of colistin (wells B4 and A4, respectively) and was therefore reported to be colistin-resistant.

intrinsically resistant to colistin (*M. morgani*, *P. mirabilis*, *P. vulgaris*, *P. stuartii*, and *S. marcescens*) and 130 isolates of various species displayed an acquired mechanism of resistance to colistin. For 87 *Klebsiella* spp. isolates, resistance to colistin was associated with various chromosomal gene changes responsible for lipopolysaccharide modifications (online Technical Appendix): 10 isolates had mutations in the PmrAB two-component system (n = 3 in *pmrA* gene, n = 7 in *pmrB* gene); 2 isolates had structural modifications in the PhoPQ two-component system (n = 1 in *phoP* gene, n = 1 in *phoQ* gene); and 75 isolates had various alterations in *mgrB* gene, the negative regulator of the PhoPQ system (online Technical Appendix). Seven nonduplicate *E. coli* isolates harbored a plasmid-mediated *mcr-1* gene. Pulsed-field gel electrophoresis revealed that isolates with identical mechanisms of colistin resistance (chromosomal or plasmid-encoded) were not clonally related (data not shown). The mechanism(s) of colistin resistance remained unknown for the 43 remaining enterobacterial isolates (online Technical Appendix). With regard to performance of the rapid polymyxin NP test with colistin-susceptible strains, the 65 colistin-susceptible isolates tested (MICs of colistin 0.12–2 µg/mL) gave negative results, except for 3 isolates (isolates FR-180, 181, and 182) for which colistin MICs were 1–2 µg/mL (just below the EUCAST breakpoint) and which gave a positive (false-positive) result (online Technical Appendix).

As expected, isolates that were intrinsically resistant to colistin (n = 5), such as *Proteus* spp., *P. stuartii* and *S. marcescens*, gave a positive test result (online Technical Appendix). Colistin-resistant enterobacterial isolates (n = 130, MICs of colistin ranging from 4 to >128 µg/mL) also gave positive results, except for 1 colistin-resistant *E. coli* isolate (isolate FR-119) for which colistin MIC was 8 µg/mL and which gave a negative (false-negative) result (online Technical Appendix).

Correlation was high between colistin resistance and positive rapid polymyxin NP test results and, conversely, colistin susceptibility and negative test results (online Technical Appendix). Sensitivity (99.3%) and specificity (95.4%) of the test were also high, compared with the standard broth microdilution method.

By reading the color change of the wells every hour, we determined that final results were obtained 2 h after incubation when the tray was incubated at 35 ± 2°C under an ambient atmosphere. However, positive results (frank color change) were obtained as early as 1 h after incubation for *Klebsiella* spp. and *E. coli* isolates. Half of the *Enterobacter* spp. isolates gave positive results within 1 h of incubation and the other half within 2 h. In addition, all resistant isolates gave positive results after 1 h of incubation when trays were incubated at 35 ± 2°C under 5% CO<sub>2</sub>. Agitating the tray did not improve the speed with which results were obtained.

The rapid polymyxin NP test results were the same whether performed with polymyxin B or with colistin (data not shown). Testing of several agar media revealed that 30% of the colistin-susceptible tested isolates gave false-positive results when bacterial colonies were recovered from acidifying media such as Drigalski, MacConkey, or bromocresol purple agar. Media that were adequate for culturing bacteria before performing the rapid polymyxin NP test were Luria Bertani agar, Mueller-Hinton agar, Columbia agar + 5% sheep blood, chocolate agar, UriSelect 4 agar, and eosin methylene blue agar.

## Discussion

The rapid polymyxin NP test is easy-to-perform, rapid, sensitive, and specific. It detects polymyxin-resistant and -susceptible isolates from any enterobacterial species, regardless of the molecular mechanism of polymyxin resistance. This test offers the possibility of detecting polymyxin resistance from bacterial cultures from infected samples or from selective media before any antimicrobial drug susceptibility testing results are obtained. Results are obtained at least 16 h sooner with this test than with the reference broth microdilution method. This test is as reliable as the reference dilution technique but much less cumbersome and is not based on diffusion of large polymyxin molecules in agar (as are the Etest and the disk-diffusion techniques), which therefore prevents false susceptibility results (15). A commercial test for research-use only is available for determining MICs of polymyxins (TREK Diagnostic Systems, Inc., Cleveland, OH, USA) (15); however, this test is adapted to testing series of isolates rather than single isolates, and results are available in 16–20 h. Sensitivity and specificity of the rapid polymyxin NP test were high (99.3% and 95.4%, respectively), making it a potential useful clinical technique.

The rapid polymyxin NP test can be performed on cultured bacteria grown on media such as Luria Bertani, Mueller-Hinton, UriSelect-4, eosin-methylene blue, blood agar, and chocolate agar. Interference may be observed with colonies grown on acidifying culture media such as Drigalski, MacConkey, and bromocresol purple agars. The test has been optimized in its present form; testing under other conditions (e.g., changes in pH indicator, inoculum size, glucose and pH indicator concentrations; preparation of the solution with non-cation-adjusted culture medium; or use of other polystyrene-containing trays) gave less optimal results.

The rapid polymyxin NP test uses commercially available polymyxin B and colistin sulfate powders, which are unspecified mixtures of chemically related compounds that differ by single amino acid changes and fatty acyl moieties. The changes of the relative proportion of the polymyxin components of these mixtures is poorly defined and a potential source of variability (21).

We believe that the rapid polymyxin NP test may be useful for first-step screening of polymyxin resistance because use of molecular-based techniques for identification of all polymyxin-resistance mechanisms cannot be foreseen. Indeed, polymyxin resistance in *Enterobacteriaceae* may be associated with many nonrelated mechanisms, some of which are known (defects in outer-membrane proteins, structural modification of lipopolysaccharide, efflux over-expression; 16) and others still unknown.

This study was subject to at least 4 limitations. First, we did not assess the ability of the rapid polymyxin NP test to detect heteroresistant isolates with low MICs for polymyxin by broth microdilution; such isolates are problematic to detect. Second, the rapid polymyxin NP test involves an orange-to-yellow color change, which is readily apparent for resistant organisms, but interpretation may require more vigilance from laboratory technicians testing organisms with low-level resistance; thus, larger scale studies in different laboratories are needed to fully evaluate the reliability of rapid polymyxin NP test performance. Third, we could not determine the mechanism for colistin resistance in 43 isolates for which several mechanisms of colistin resistance may exist and may be expressed in various ways. Last, we did not evaluate the rapid polymyxin NP test in species of bacteria with metabolic pathways that differed from those of *Enterobacteriaceae*; further work is needed to adapt the rapid polymyxin NP test to detection of polymyxin-resistant *Pseudomonas aeruginosa* and *Acinetobacter baumannii*, which have different metabolic pathways.

The rapid polymyxin NP test can be used to determine susceptibility or resistance to polymyxins in countries facing endemic spread of carbapenemase producers and for which polymyxins are last-resort drugs (22). The test can also rapidly identify carriers of polymyxin-resistant isolates, leading to rapid implementation of adequate hygiene measures to control their spread. This test can also support the development of novel polymyxin-like molecules, facilitating patient enrollment in pivotal clinical trials.

### Acknowledgments

We thank the following colleagues, who provided polymyxin-resistant isolates: Adrian Brink, Jean-Winoc Decousser, Javier Fernandez-Dominguez, Jean-Yves Madec, Maria-Virginia Villegas, and Mesut Yilmaz.

This work was financed by the University of Fribourg, Switzerland. An international patent form, corresponding to this test, has been filed on behalf of the University of Fribourg.

Dr. Nordmann is professor of medicine, chair of the Microbiology Unit, Department of Medicine, Faculty of Science, University of Fribourg, Switzerland. His main field of research interest includes the genetics, biochemistry, and molecular epidemiology

of resistance in gram-negative bacteria, and the development of rapid diagnostic tests for detection of emerging antibiotic resistance traits.

### References

- Nordmann P, Naas T, Poirel L. Global spread of carbapenemase-producing *Enterobacteriaceae*. *Emerg Infect Dis*. 2011;17:1791–8. <http://dx.doi.org/10.3201/eid1710.110655>
- Falagas ME, Kasiakou SK. Colistin: the revival of polymyxins for the management of multidrug-resistant gram-negative bacterial infections. *Clin Infect Dis*. 2005;40:1333–41. <http://dx.doi.org/10.1086/429323>
- Monaco M, Giani T, Raffone M, Arena F, Garcia-Fernandez A, Pollini S, et al. Colistin resistance superimposed to endemic carbapenem-resistant *Klebsiella pneumoniae*: a rapidly evolving problem in Italy, November 2013 to April 2014. *Eur Surveill*. 2014;19:pii:20939.
- Hindler JA, Humphries RM. Colistin MIC variability by method for contemporary clinical isolates of multidrug-resistant Gram-negative bacilli. *J Clin Microbiol*. 2013;51:1678–84. <http://dx.doi.org/10.1128/JCM.03385-12>
- Tan TY, Ng SY. Comparison of Etest, Vitek and agar dilution for susceptibility testing of colistin. *Clin Microbiol Infect*. 2007;13:541–4. <http://dx.doi.org/10.1111/j.1469-0691.2007.01708.x>
- Olaitan AO, Morand S, Rolain JM. Mechanisms of polymyxin resistance: acquired and intrinsic resistance in bacteria. *Front Microbiol*. 2014;5:643. <http://dx.doi.org/10.3389/fmicb.2014.00643>
- Liu YY, Wang Y, Walsh TR, Yi LX, Zhang R, Spencer J, et al. Emergence of plasmid-mediated colistin resistance mechanism MCR-1 in animals and human beings in China: a microbiological and molecular biological study. *Lancet Infect Dis*. 2016;16:161–8. [http://dx.doi.org/10.1016/S1473-3099\(15\)00424-7](http://dx.doi.org/10.1016/S1473-3099(15)00424-7)
- Arcilla MS, van Hattem JM, Matamoros S, Melles DC, Penders J, de Jong MD, et al. Dissemination of the *mcr-1* colistin resistance gene. *Lancet Infect Dis*. 2016;16:147–9. [http://dx.doi.org/10.1016/S1473-3099\(15\)00541-1](http://dx.doi.org/10.1016/S1473-3099(15)00541-1)
- Malhotra-Kumar S, Xavier BB, Das AJ, Lammens C, Hoang HT, Pham TP, et al. Colistin resistance gene *mcr-1* harboured on a multidrug resistant plasmid. *Lancet Infect Dis*. 2016; 2016;16:283–4. [http://dx.doi.org/10.1016/S1473-3099\(16\)00012-8](http://dx.doi.org/10.1016/S1473-3099(16)00012-8)
- Poirel L, Kieffer N, Liassine N, Thanh D, Nordmann P. Plasmid-mediated carbapenem and colistin resistance in a clinical isolate of *Escherichia coli*. *Lancet Infect Dis*. 2016;16:281. [http://dx.doi.org/10.1016/S1473-3099\(16\)00006-2](http://dx.doi.org/10.1016/S1473-3099(16)00006-2)
- Hugh R, Leifson E. The taxonomic significance of fermentative versus oxidative metabolism of carbohydrates by various gram negative bacteria. *J Bacteriol*. 1953;66:24–6.
- Jayol A, Nordmann P, Brink A, Poirel L. Heteroresistance to colistin in *Klebsiella pneumoniae* associated with alterations in the PhoPQ regulatory system. *Antimicrob Agents Chemother*. 2015;59:2780–4. <http://dx.doi.org/10.1128/AAC.05055-14>
- Clinical and Laboratory Standards Institute. Methods for dilution of antimicrobial susceptibility tests for bacteria that grow aerobically: approved standard. 9th ed. Document M07–A9. Wayne (PA): The Institute; 2012.
- Clinical and Laboratory Standards Institute. Performance standards for antimicrobial susceptibility testing: 24th informational supplement. Document M100–S24. Wayne (PA): The Institute; 2014.
- Humphries RM. Susceptibility testing of the polymyxins: where are we now? *Pharmacotherapy*. 2015;35:22–7. <http://dx.doi.org/10.1002/phar.1505>
- European Committee on Antimicrobial Susceptibility Testing. Breakpoints tables for interpretation of MICs and zone diameters.

Version 2.0. 2014 [cited 2016 Apr 5]. [http://www.eucast.org/fileadmin/src/media/PDFs/EUCAST\\_files/Breakpoint\\_tables/v\\_6.0\\_Breakpoint\\_table.pdf](http://www.eucast.org/fileadmin/src/media/PDFs/EUCAST_files/Breakpoint_tables/v_6.0_Breakpoint_table.pdf)


17. Jayol A, Poirel L, Brink A, Villegas MV, Yilmaz M, Nordmann P. Resistance to colistin associated with a single amino acid change in protein PmrB among *Klebsiella pneumoniae* isolates of worldwide origin. *Antimicrob Agents Chemother*. 2014;58:4762–6. <http://dx.doi.org/10.1128/AAC.00084-14>
18. Jayol A, Poirel L, Villegas MV, Nordmann P. Modulation of *mgrB* gene expression as a source of colistin resistance in *Klebsiella oxytoca*. *Int J Antimicrob Agents*. 2015;46:108–10. <http://dx.doi.org/10.1016/j.ijantimicag.2015.02.015>
19. Poirel L, Jayol A, Bontron S, Villegas MV, Ozdamar M, Turkoglu S, et al. The *mgrB* gene as a key target for acquired resistance to colistin in *Klebsiella pneumoniae*. *J Antimicrob Chemother*. 2015;70:75–80. <http://dx.doi.org/10.1093/jac/dku323>
20. Carrër A, Lassel L, Fortineau N, Mansouri M, Anguel N, Richard C, et al. Outbreak of CTX-M-15-producing *Klebsiella pneumoniae* in the intensive care unit of a French hospital. *Microb Drug Resist*. 2009;15:47–54. <http://dx.doi.org/10.1089/mdr.2009.0868>
21. He J, Ledesma KR, Lam WY, Figueroa DA, Lim TP, Chow DS, et al. Variability of polymyxin B major components in commercial formulations. *Int J Antimicrob Agents*. 2010;35:308–10. <http://dx.doi.org/10.1016/j.ijantimicag.2009.11.005>
22. Capone A, Giannella M, Fortini D, Giordano A, Meledandri M, Ballardini M, et al. High rate of colistin resistance among patients with carbapenem-resistant *Klebsiella pneumoniae* infection accounts for an excess of mortality. *Clin Microbiol Infect*. 2013;19:E23–30. <http://dx.doi.org/10.1111/1469-0691.12070>

Address for correspondence: Patrice Nordmann, Medical and Molecular Microbiology, Department of Medicine, Faculty of Science, University of Fribourg, Rue Albert Gockel 3, CH-1700 Fribourg, Switzerland; email: [patrice.nordmann@unifr.ch](mailto:patrice.nordmann@unifr.ch)

## December 2015: Zoonotic Infections

### Including:

- Identifying and Reducing Remaining Stocks of Rinderpest Virus
- Opportunistic Pulmonary *Bordetella hinzii* Infection after Avian Exposure
- Zoonotic Leprosy in the Southeastern United States
- Infection Risk for Persons Exposed to Highly Pathogenic Avian Influenza A H5 Virus–Infected Birds, United States, December 2014–March 2015
- High Prevalence of Intermediate *Leptospira* spp. DNA in Febrile Humans From Urban and Rural Ecuador
- Methicillin-Resistant *Staphylococcus aureus* Prevalence among Captive Chimpanzees, Texas, USA, 2012

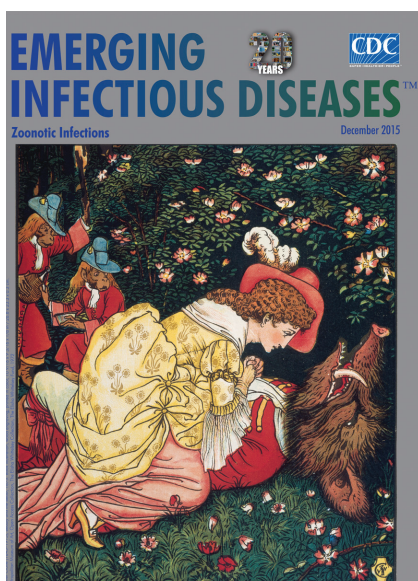

- Influenza A(H6N1) Virus in Dogs, Taiwan
- Novel *Waddlia* Intracellular Bacterium in *Artibeus intermedius* Fruit Bats, Mexico
- Tembusu-Related Flavivirus in Ducks, Thailand
- Japanese Macaques (*Macaca fuscata*) as Natural Reservoir of *Bartonella quintana*
- Increased Number of Human Cases of Influenza Virus A(H5N1) Infection, Egypt, 2014–15
- Replication Capacity of Avian Influenza A(H9N2) Virus in Pet Birds, Chickens, and Mammals, Bangladesh
- Pyrethroid and DDT Resistance and Organophosphate Susceptibility among *Anopheles* spp. Mosquitoes, Western Kenya

<http://wwwnc.cdc.gov/eid/articles/issue/21/12/table-of-contents>

# Rapid Detection of Polymyxin Resistance in *Enterobacteriaceae*

## Technical Appendix

**Technical Appendix Table.** Rapid polymyxin NP test results for polymyxin-resistant isolates with intrinsic resistance, chromosome- and plasmid-mediated acquired resistance, and for polymyxin-susceptible isolates\*

Strain	Species	Origin	Polymyxin phenotype	Rapid polymyxin NP test	Colistin MIC, µg/mL	Resistance mechanism	Strain reference
Strains intrinsically resistant to colistin							
FR-01†	<i>M. organii</i>	France	R	+	>128	NA	Positive control
FR-02	<i>P. mirabilis</i>	France	R	+	>128	NA	Unpublished
FR-03	<i>P. vulgaris</i>	France	R	+	>128	NA	Unpublished
FR-04	<i>P. stuartii</i>	France	R	+	>128	NA	Unpublished
FR-05	<i>S. marcescens</i>	France	R	+	>128	NA	Unpublished
Strains resistant to colistin with an identified mechanism of resistance							
FR-06	<i>K. pneumoniae</i>	Colombia	R	+	32	PmrA G53C	Unpublished
FR-07	<i>K. pneumoniae</i>	France	R	+	32	PmrA G53S	Unpublished
FR-08	<i>K. pneumoniae</i>	South Africa	R	+	128	PmrA G53S	Unpublished
FR-09	<i>K. pneumoniae</i>	Turkey	R	+	32	PmrB L17Q	Unpublished
FR-10	<i>K. pneumoniae</i>	South Africa	R	+	16	PmrB T157P	Isolate AF1b (1)
FR-11	<i>K. pneumoniae</i>	Colombia	R	+	32	PmrB T157P	Isolate C3 (1)
FR-12	<i>K. pneumoniae</i>	Colombia	R	+	16	PmrB T157P	Isolate C19 (1)
FR-13	<i>K. pneumoniae</i>	Turkey	R	+	32	PmrB T157P	Isolate T2 (1)
FR-14	<i>K. pneumoniae</i>	Turkey	R	+	32	PmrB T157P	Unpublished
FR-15	<i>K. pneumoniae</i>	Turkey	R	+	16	PmrB T157P	Unpublished
FR-16	<i>K. pneumoniae</i>	South Africa	HR	+	128	PhoP D191Y±deletion of 25 nt	Isolate Kp75 (2)
FR-17	<i>K. pneumoniae</i>	Turkey	R	+	>128	PhoQ R16C	Unpublished
FR-18	<i>K. pneumoniae</i>	Turkey	R	+	32	MgrB W20R	Unpublished
FR-19	<i>K. pneumoniae</i>	France	R	+	32	MgrB M27K	Unpublished
FR-20	<i>K. pneumoniae</i>	France	R	+	64	MgrB C39Y	Unpublished
FR-21	<i>K. pneumoniae</i>	France	R	+	32	MgrB N42Y and K43I	Unpublished
FR-22	<i>K. pneumoniae</i>	Angola	R	+	64	MgrB I45T	Unpublished
FR-23	<i>K. pneumoniae</i>	France	R	+	64	MgrB P46S	Unpublished
FR-24	<i>K. pneumoniae</i>	Turkey	R	+	4	MgrB W47R	Unpublished
FR-25	<i>K. pneumoniae</i>	France	R	+	128	MgrB truncated (2 aa)	Unpublished
FR-26	<i>K. pneumoniae</i>	France	R	+	128	MgrB truncated (2 aa)	Unpublished
FR-27	<i>K. pneumoniae</i>	France	R	+	64	MgrB truncated (2 aa)	Unpublished
FR-28	<i>K. pneumoniae</i>	France	R	+	128	MgrB truncated (27 aa)	Unpublished
FR-29	<i>K. pneumoniae</i>	France	R	+	32	MgrB truncated (27 aa)	Unpublished
FR-30	<i>K. pneumoniae</i>	France	R	+	>128	MgrB truncated (27 aa)	Isolate Sa (3)
FR-31	<i>K. pneumoniae</i>	France	R	+	64	MgrB truncated (29 aa)	Isolate 15I5 (3)
FR-32	<i>K. pneumoniae</i>	Switzerland	R	+	64	MgrB truncated (29 aa)	Unpublished
FR-33	<i>K. pneumoniae</i>	Spain	R	+	32	MgrB truncated (29 aa)	Unpublished
FR-34	<i>K. pneumoniae</i>	South Africa	R	+	64	MgrB truncated (29 aa)	Unpublished
FR-35	<i>K. pneumoniae</i>	South Africa	R	+	128	MgrB truncated (29 aa)	Unpublished
FR-36	<i>K. pneumoniae</i>	Colombia	R	+	128	MgrB truncated (29 aa)	Isolate C11 (3)
FR-37	<i>K. pneumoniae</i>	France	R	+	32	MgrB truncated (32 aa)	Unpublished
FR-38	<i>K. pneumoniae</i>	Turkey	R	+	32	MgrB truncated (46 aa)	Unpublished
FR-39	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> IS 1R between +21 and +22	Unpublished
FR-40	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> ISEcp1/ <i>bla</i> CTX-M-15 between +21 and +22	Unpublished
FR-41	<i>K. pneumoniae</i>	France	R	+	32	<i>mgrB</i> IS 102-like between +36 and +37	Unpublished

Strain	Species	Origin	Polymyxin phenotype	Rapid polymyxin NP test	Colistin MIC, µg/mL	Resistance mechanism	Strain reference
FR-42	<i>K. pneumoniae</i>	France	R	+	>128	<i>mgrB</i> IS102-like between +44 and +45	Unpublished
FR-43	<i>K. pneumoniae</i>	South Africa	R	+	64	<i>mgrB</i> IS903b between +44 and +45	Unpublished
FR-44	<i>K. pneumoniae</i>	Turkey	R	+	64	<i>mgrB</i> IS2 between +44 and +45	Unpublished
FR-45	<i>K. pneumoniae</i>	Turkey	R	+	128	<i>mgrB</i> IS1R between +44 and +45	Unpublished
FR-46	<i>K. pneumoniae</i>	South Africa	R	+	64	<i>mgrB</i> IS1R between +61 and +62	Unpublished
FR-47	<i>K. pneumoniae</i>	Turkey	R	+	64	<i>mgrB</i> IS903b-like between +69 and +70	Unpublished
FR-48	<i>K. pneumoniae</i>	Spain	R	+	128	<i>mgrB</i> IS903-like between +70 and +71	Unpublished
FR-49	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> IS5-like between +74 and 75	Isolate 1118 (3)
FR-50	<i>K. pneumoniae</i>	France	R	+	16	<i>mgrB</i> IS5-like between +74 and 75	Isolate 20C9 (3)
FR-51	<i>K. pneumoniae</i>	Colombia	R	+	64	<i>mgrB</i> IS5-like between +74 and 75	Isolate C9 (3)
FR-52	<i>K. pneumoniae</i>	Turkey	R	+	128	<i>mgrB</i> IS5-like between +74 and 75	Isolate T1b (3)
FR-53	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> IS5-like between +74 and +75	Unpublished
FR-54	<i>K. pneumoniae</i>	Colombia	R	+	128	<i>mgrB</i> ISKpn13 between +74 and +75	Isolate C21 (3)
FR-55	<i>K. pneumoniae</i>	France	R	+	>128	<i>mgrB</i> ISKpn26-like between +74 and +75	Unpublished
FR-56	<i>K. pneumoniae</i>	Spain	R	+	64	<i>mgrB</i> ISKpn26-like between +74 and +75	Unpublished
FR-57	<i>K. pneumoniae</i>	Turkey	R	+	128	<i>mgrB</i> ISKpn26-like between +74 and +75	Unpublished
FR-58	<i>K. pneumoniae</i>	South Africa	R	+	32	<i>mgrB</i> ISKpn26-like between +76 and +77	Unpublished
FR-59	<i>K. pneumoniae</i>	South Africa	R	+	32	<i>mgrB</i> IS903B between +76 and +77	Unpublished
FR-60	<i>K. pneumoniae</i>	Turkey	R	+	128	<i>mgrB</i> ISKpn14 between +77 and +78	Unpublished
FR-61	<i>K. pneumoniae</i>	France	R	+	128	<i>mgrB</i> duplication 19 nt between +84 and +85	Unpublished
FR-62	<i>K. pneumoniae</i>	South Africa	R	+	64	<i>mgrB</i> IS903b between +94 and +95	Unpublished
FR-63	<i>K. pneumoniae</i>	France	R	+	128	<i>mgrB</i> IS1R-like between +116 and +117	Unpublished
FR-64	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> IS903b-like between +116 and +117	Unpublished
FR-65	<i>K. pneumoniae</i>	France	R	+	8	<i>mgrB</i> IS1R between +118 and +119	Unpublished
FR-66	<i>K. pneumoniae</i>	Turkey	R	+	32	<i>mgrB</i> IS1R between +123 and +124	Unpublished
FR-67	<i>K. pneumoniae</i>	France	R	+	128	<i>mgrB</i> ISKpn26-like between +125 and +126	Unpublished
FR-68	<i>K. pneumoniae</i>	Colombia	R	+	64	<i>mgrB</i> ISKpn14 between +127 and +128	Isolate C22 (3)
FR-69	<i>K. pneumoniae</i>	South Africa	R	+	128	<i>mgrB</i> IS1R between +131 and +132	Unpublished
FR-70	<i>K. pneumoniae</i>	Colombia	R	+	128	<i>mgrB</i> IS10R between -26 and -27	Isolate C1 (3)
FR-71	<i>K. pneumoniae</i>	Turkey	R	+	64	<i>mgrB</i> ISKpn14 between -27 and -28	Unpublished
FR-72	<i>K. pneumoniae</i>	Turkey	R	+	64	<i>mgrB</i> ISKpn14 between -28 and -29	Unpublished
FR-73	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> IS1R between -36 and -37	Unpublished
FR-74	<i>K. pneumoniae</i>	France	R	+	128	<i>mgrB</i> IS1R between -45 and -46	Unpublished



Strain	Species	Origin	Polymyxin phenotype	Rapid polymyxin NP test	Colistin MIC, µg/mL	Resistance mechanism	Strain reference
FR-75	<i>K. pneumoniae</i>	France	R	+	64	<i>mgrB</i> IS1R between -45 and -46	Unpublished
FR-76	<i>K. pneumoniae</i>	Turkey	R	+	32	<i>mgrB</i> IS1R between -45 and -46	Unpublished
FR-77	<i>K. pneumoniae</i>	France	R	+	32	<i>mgrB</i> ISKpn14-like between -45 and -46	Unpublished
FR-78	<i>K. pneumoniae</i>	France	R	+	32	<i>mgrB</i> IS1R between -61 and -62	Unpublished
FR-79	<i>K. pneumoniae</i>	France	R	+	16	Full <i>mgrB</i> gene deletion	Unpublished
FR-80	<i>K. pneumoniae</i>	France	R	+	>128	Full <i>mgrB</i> gene deletion	Unpublished
FR-81	<i>K. pneumoniae</i>	South Africa	R	+	64	Full <i>mgrB</i> gene deletion	Unpublished
FR-82	<i>K. pneumoniae</i>	South Africa	R	+	64	Full <i>mgrB</i> gene deletion	Unpublished
FR-83	<i>K. pneumoniae</i>	France	R	+	32	Deletion nt 23 <i>mgrB</i>	Unpublished
FR-84	<i>K. pneumoniae</i>	South Africa	R	+	128	Deletion nt 70 <i>mgrB</i> and substitution nt 73	Unpublished
FR-85	<i>K. pneumoniae</i>	France	R	+	64	Deletion nt 74 <i>mgrB</i>	Unpublished
FR-86	<i>K. pneumoniae</i>	Spain	R	+	64	Deletion nt 100 <i>mgrB</i>	Unpublished
FR-87	<i>K. pneumoniae</i>	Spain	R	+	16	Deletion nt 100 <i>mgrB</i>	Unpublished
FR-88	<i>K. pneumoniae</i>	Turkey	R	+	>128	Deletion nt 22 to 32 <i>mgrB</i>	Unpublished
FR-89	<i>K. pneumoniae</i>	Colombia	R	+	>128	Deletion nt 23 to 33 <i>mgrB</i>	Unpublished
FR-90	<i>K. pneumoniae</i>	South Africa	R	+	64	Deletion nt 30 et 31 <i>mgrB</i>	Unpublished
FR-91	<i>K. pneumoniae</i>	South Africa	R	+	64	Deletion nt 48 to 57 <i>mgrB</i>	Unpublished
FR-92	<i>K. oxytoca</i>	Colombia	R	+	64	MgrB ISKpn26-like between -38 and -39	Isolate C24 (4)
FR-93	<i>E. coli</i>	Switzerland	R	+	4	Plasmid-mediated <i>mcr-1</i> gene	Isolate KRI (5)
FR-94	<i>E. coli</i>	South Africa	R	+	16	Plasmid-mediated <i>mcr-1</i> gene	Unpublished
FR-95	<i>E. coli</i>	South Africa	R	+	16	Plasmid-mediated <i>mcr-1</i> gene	Unpublished
FR-96	<i>E. coli</i>	South Africa	R	+	16	Plasmid-mediated <i>mcr-1</i> gene	Unpublished
FR-97	<i>E. coli</i>	South Africa	R	+	16	Plasmid-mediated <i>mcr-1</i> gene	Unpublished
FR-98	<i>E. coli</i>	South Africa	R	+	8	Plasmid-mediated <i>mcr-1</i> gene	Unpublished
FR-99	<i>E. coli</i>	South Africa	R	+	16	Plasmid-mediated <i>mcr-1</i> gene	Unpublished
Strains resistant to colistin without an identified mechanism of resistance							
FR-100	<i>K. pneumoniae</i>	France	R	+	32	Unknown	Unpublished
FR-101	<i>K. pneumoniae</i>	France	R	+	32	Unknown	Unpublished
FR-102	<i>K. pneumoniae</i>	France	R	+	32	Unknown	Unpublished
FR-103	<i>K. pneumoniae</i>	France	R	+	4	Unknown	Unpublished
FR-104	<i>K. pneumoniae</i>	France	R	+	64	Unknown	Unpublished
FR-105	<i>K. pneumoniae</i>	France	R	+	128	Unknown	Unpublished
FR-106	<i>K. pneumoniae</i>	France	R	+	>128	Unknown	Unpublished
FR-107	<i>K. pneumoniae</i>	Spain	R	+	32	Unknown	Unpublished
FR-108	<i>K. pneumoniae</i>	Spain	R	+	32	Unknown	Unpublished
FR-109	<i>K. pneumoniae</i>	South Africa	R	+	64	Unknown	Unpublished
FR-110	<i>K. pneumoniae</i>	South Africa	R	+	>128	Unknown	Unpublished
FR-111	<i>K. pneumoniae</i>	South Africa	R	+	32	Unknown	Unpublished
FR-112	<i>K. pneumoniae</i>	Colombia	R	+	>128	Unknown	Unpublished
FR-113	<i>K. pneumoniae</i>	Colombia	R	+	>128	Unknown	Unpublished
FR-114	<i>K. pneumoniae</i>	Colombia	R	+	64	Unknown	Unpublished
FR-115	<i>K. pneumoniae</i>	Colombia	R	+	64	Unknown	Unpublished
FR-116	<i>K. pneumoniae</i>	Colombia	R	+	32	Unknown	Unpublished
FR-117	<i>K. pneumoniae</i>	Turkey	R	+	32	Unknown	Unpublished
FR-118	<i>K. pneumoniae</i>	Turkey	R	+	4	Unknown	Unpublished
FR-119	<i>E. coli</i>	France	R	-	8	Unknown	Unpublished
FR-120	<i>E. coli</i>	France	R	+	8	Unknown	Unpublished
FR-121	<i>E. coli</i>	France	R	+	4	Unknown	Unpublished
FR-122	<i>E. cloacae</i>	Colombia	R	+	32	Unknown	Unpublished

Strain	Species	Origin	Polymyxin phenotype	Rapid polymyxin NP test	Colistin MIC, µg/mL	Resistance mechanism	Strain reference
FR-123	<i>E. cloacae</i>	Colombia	R	+	>128	Unknown	Unpublished
FR-124	<i>E. cloacae</i>	Colombia	R	+	>128	Unknown	Unpublished
FR-125	<i>E. cloacae</i>	France	R	+	>128	Unknown	Unpublished
FR-126	<i>E. cloacae</i>	France	R	+	>128	Unknown	Unpublished
FR-127	<i>E. cloacae</i>	France	R	+	64	Unknown	Unpublished
FR-128	<i>E. cloacae</i>	France	R	+	64	Unknown	Unpublished
FR-129	<i>E. cloacae</i>	France	R	+	32	Unknown	Unpublished
FR-130	<i>E. cloacae</i>	France	R	+	>128	Unknown	Unpublished
FR-131	<i>E. cloacae</i>	France	R	+	16	Unknown	Unpublished
FR-132	<i>E. cloacae</i>	France	R	+	>128	Unknown	Unpublished
FR-133	<i>E. cloacae</i>	France	R	+	>128	Unknown	Unpublished
FR-134	<i>E. asburiae</i>	France	R	+	>128	Unknown	Unpublished
FR-135	<i>H. alvei</i>	France	R	+	>128	Unknown	Unpublished
Strains susceptible to colistin							
FR-136 <sup>†</sup>	<i>E. coli</i>	USA	S	-	0.25	NA	ATCC 25922, Negative control
FR-137	<i>E. coli</i>	Colombia	S	-	0.12	NA	Unpublished
FR-138	<i>E. coli</i>	Colombia	S	-	0.12	NA	Unpublished
FR-139	<i>E. coli</i>	Switzerland	S	-	0.12	NA	Unpublished
FR-140	<i>E. coli</i>	Switzerland	S	-	0.12	NA	Unpublished
FR-141	<i>E. coli</i>	Switzerland	S	-	0.12	NA	Unpublished
FR-142	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-143	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-144	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-145	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-146	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-147	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-148	<i>E. coli</i>	France	S	-	0.12	NA	Unpublished
FR-149	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-150	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-151	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-152	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-153	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-154	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-155	<i>E. coli</i>	France	S	-	0.25	NA	Unpublished
FR-156	<i>K. pneumoniae</i>	USA	S	-	0.12	NA	ATCC 53153
FR-157	<i>K. pneumoniae</i>	South Africa	S	-	0.25	NA	Isolate AF1a (1)
FR-158	<i>K. pneumoniae</i>	Turkey	S	-	0.12	NA	Isolate T1a (3)
FR-159	<i>K. pneumoniae</i>	Spain	S	-	0.12	NA	Unpublished
FR-160	<i>K. pneumoniae</i>	Spain	S	-	0.25	NA	Unpublished
FR-161	<i>K. pneumoniae</i>	Spain	S	-	0.25	NA	Unpublished
FR-162	<i>K. pneumoniae</i>	Spain	S	-	0.5	NA	Unpublished
FR-163	<i>K. pneumoniae</i>	Spain	S	-	0.5	NA	Unpublished
FR-164	<i>K. pneumoniae</i>	Spain	S	-	0.5	NA	Unpublished
FR-165	<i>K. pneumoniae</i>	Spain	S	-	1	NA	Unpublished
FR-166	<i>K. pneumoniae</i>	Colombia	S	-	0.12	NA	Unpublished
FR-167	<i>K. pneumoniae</i>	Colombia	S	-	0.12	NA	Unpublished
FR-168	<i>K. pneumoniae</i>	Colombia	S	-	0.12	NA	Unpublished
FR-169	<i>K. pneumoniae</i>	Colombia	S	-	0.25	NA	Unpublished
FR-170	<i>K. pneumoniae</i>	Colombia	S	-	0.5	NA	Unpublished
FR-171	<i>K. pneumoniae</i>	Switzerland	S	-	0.12	NA	Unpublished
FR-172	<i>K. pneumoniae</i>	Switzerland	S	-	0.12	NA	Unpublished
FR-173	<i>K. pneumoniae</i>	Switzerland	S	-	0.25	NA	Unpublished
FR-174	<i>K. pneumoniae</i>	France	S	-	0.12	NA	Unpublished
FR-175	<i>K. pneumoniae</i>	France	S	-	0.12	NA	Unpublished
FR-176	<i>K. pneumoniae</i>	France	S	-	0.12	NA	Unpublished
FR-177	<i>K. pneumoniae</i>	France	S	-	0.25	NA	Unpublished
FR-178	<i>K. pneumoniae</i>	France	S	-	0.25	NA	Unpublished
FR-179	<i>K. pneumoniae</i>	France	S	-	0.25	NA	Unpublished
FR-180	<i>K. pneumoniae</i>	France	S	+	1	NA	Unpublished
FR-181	<i>K. pneumoniae</i>	France	S	+	2	NA	Unpublished
FR-182	<i>K. pneumoniae</i>	France	S	+	2	NA	Unpublished
FR-183	<i>K. oxytoca</i>	France	S	-	0.12	NA	Unpublished
FR-184	<i>K. oxytoca</i>	France	S	-	0.12	NA	Unpublished
FR-185	<i>K. oxytoca</i>	France	S	-	0.25	NA	Unpublished

Strain	Species	Origin	Polymyxin phenotype	Rapid polymyxin NP test	Colistin MIC, µg/mL	Resistance mechanism	Strain reference
FR-186	<i>E. cloacae</i>	Colombia	S	–	0.12	NA	Unpublished
FR-187	<i>E. cloacae</i>	Colombia	S	–	0.12	NA	Unpublished
FR-188	<i>E. cloacae</i>	France	S	–	0.12	NA	Unpublished
FR-189	<i>E. cloacae</i>	France	S	–	0.12	NA	Unpublished
FR-190	<i>E. cloacae</i>	France	S	–	0.25	NA	Unpublished
FR-191	<i>E. cloacae</i>	France	S	–	0.5	NA	Unpublished
FR-192	<i>E. aerogenes</i>	Switzerland	S	–	0.12	NA	Unpublished
FR-193	<i>E. aerogenes</i>	France	S	–	0.12	NA	Unpublished
FR-194	<i>E. aerogenes</i>	France	S	–	0.25	NA	Unpublished
FR-195	<i>C. freundii</i>	France	S	–	0.12	NA	Unpublished
FR-196	<i>C. freundii</i>	France	S	–	0.25	NA	Unpublished
FR-197	<i>C. freundii</i>	Colombia	S	–	0.25	NA	Unpublished
FR-198	<i>C. koseri</i>	France	S	–	0.12	NA	Unpublished
FR-199	<i>C. koseri</i>	France	S	–	0.12	NA	Unpublished
FR-200	<i>C. koseri</i>	France	S	–	0.25	NA	Unpublished

\*HR, heteroresistant; NA, not applicable; R, resistant; S, susceptible.

†Strains FR-01 and FR-136 were used as positive and negative controls, respectively.

## References

- Jayol A, Poirel L, Brink A, Villegas MV, Yilmaz M, Nordmann P. Resistance to colistin associated with a single amino acid change in protein PmrB among *Klebsiella pneumoniae* isolates of worldwide origin. *Antimicrob Agents Chemother.* 2014;58:4762–6. [PubMed](#)  
<http://dx.doi.org/10.1128/AAC.00084-14>
- Jayol A, Nordmann P, Brink A, Poirel L. Heteroresistance to colistin in *Klebsiella pneumoniae* associated with alterations in the PhoPQ regulatory system. *Antimicrob Agents Chemother.* 2015;59:2780–4. [PubMed](#) <http://dx.doi.org/10.1128/AAC.05055-14>
- Poirel L, Jayol A, Bontron S, Villegas MV, Ozdamar M, Turkoglu S, et al. The *mgrB* gene as a key target for acquired resistance to colistin in *Klebsiella pneumoniae*. *J Antimicrob Chemother.* 2015;70:75–80. [PubMed](#) <http://dx.doi.org/10.1093/jac/dku323>
- Jayol A, Poirel L, Villegas MV, Nordmann P. Modulation of *mgrB* gene expression as a source of colistin resistance in *Klebsiella oxytoca*. *Int J Antimicrob Agents.* 2015;46:108–10. [PubMed](#)  
<http://dx.doi.org/10.1016/j.ijantimicag.2015.02.015>
- Poirel L, Kieffer N, Liassine N, Thanh D, Nordmann P. Plasmid-mediated carbapenem and colistin resistance in a clinical isolate of *Escherichia coli*. *Lancet Infect Dis.* 2016;16:281. [PubMed](#)  
[http://dx.doi.org/10.1016/S1473-3099\(16\)00006-2](http://dx.doi.org/10.1016/S1473-3099(16)00006-2)